



1. Experimental Procedure

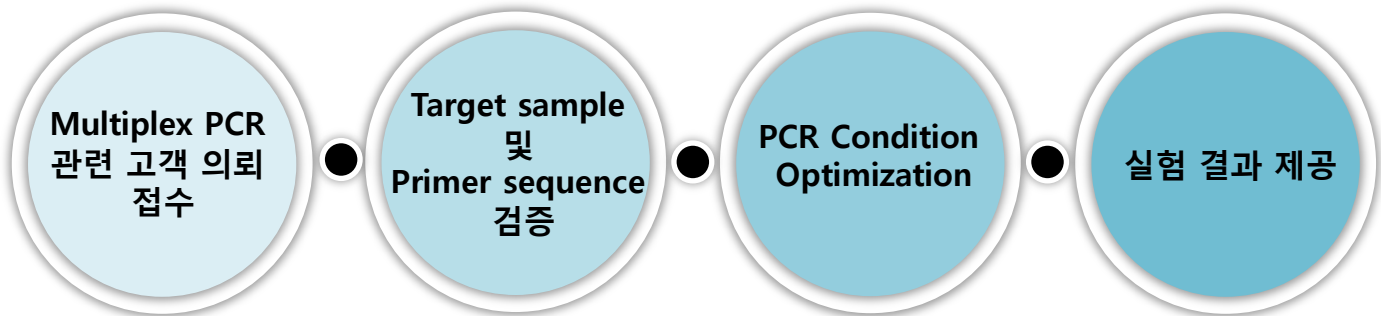

2. Experimental Results

Sample type	Protozoa Penal
Sample amount	10 ³ ~ 10 ² Copy
PCR condition	Pre-denaturation : 95 °C (3min) Denaturation : 95 °C (30sec) Annealing : 60 °C (30sec) Extension : 72 °C (30sec) to 45cycles

Multiplex PCR PreMix

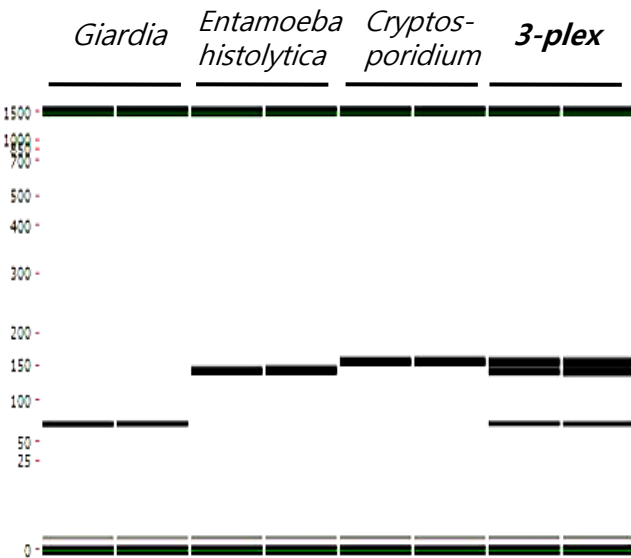

Supplier A, PCR Master Mix

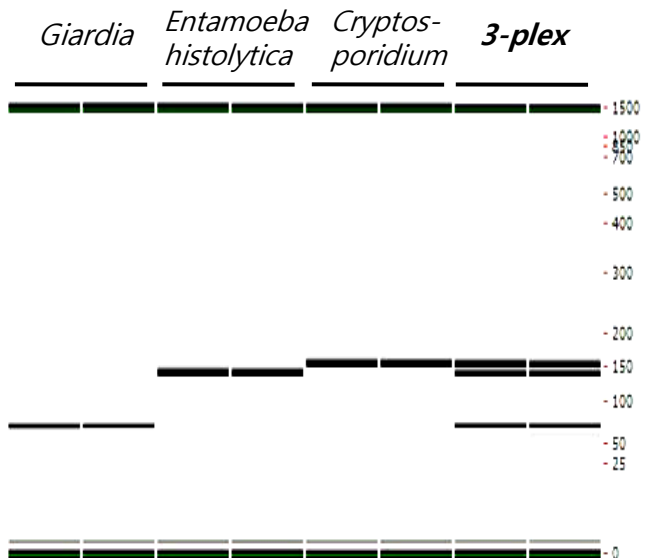


Fig 1. Comparison of amplification quality using Lapchip[®] between AccuPower Multiplex PCR PreMix and other supplier's master mixture.